

A Protocol Guide for the *N. crassa* Yeast Artificial Chromosome Library

M. B. Centola(1), *C. T. Yamashiro*(2), *L. S. Martel*(3), *J. C. Royer*(4) and *T. J. Schmidhauser*(5).
 (1)*XOMA Inc., 1545 17th St., Santa Monica, CA 90404-3401.* (2)*Department of Biological Sciences, Stanford University, Stanford, CA 94305-5020.* (3)*Department of Microbiology and Molecular Genetics, University of California, Los Angeles, Los Angeles, CA 90024-1570.* (4)*Novo Nordisk Biotech, Inc., 1445 Drew Ave., Davis, CA 95616-4880.* (5)*Department of Medical Biochemistry, Southern Illinois University, Carbondale, IL 62901-4413.*

A yeast artificial chromosome (YAC) library of *Neurospora crassa* strain 74-OR23-1A has been constructed. This library has been used to clone 750 kb of contiguous DNA sequences from the centromere region of linkage group VII (M. Centola and J. Carbon. 1994. *Mol. Cell. Biol.* 14:1510-1519). The purpose of this article is explicitly to outline procedures that have been developed for library screening and chromosome walking. The library was constructed in the YAC vector pYAC4 (Burke et al. 1987. *Science* 236:806-812). This vector contains *Tetrahymena* telomere, *Saccharomyces cerevisiae* CEN4, and ARS1 DNA sequences that specify in cis full telomeric, centromeric and replication functions in *S. cerevisiae*. YAC clones are maintained in approximately single copy as highly stable linear "artificial chromosomes". pYAC4 also encodes the yeast TRP1 and URA3 for selection of YAC clones in the yeast host strain AB1380 MATa ade2-1 can1-100 lys2-1 trp1 ura3 his5 {psi+} (Burke et al. 1987). A model of a YAC clone is shown in Fig. 1. The *N. crassa* YAC library contains 2204 clones, with an average insert size of 170 kb. Inserts from 40 clones have been characterized and range in size from 75-260 kb.

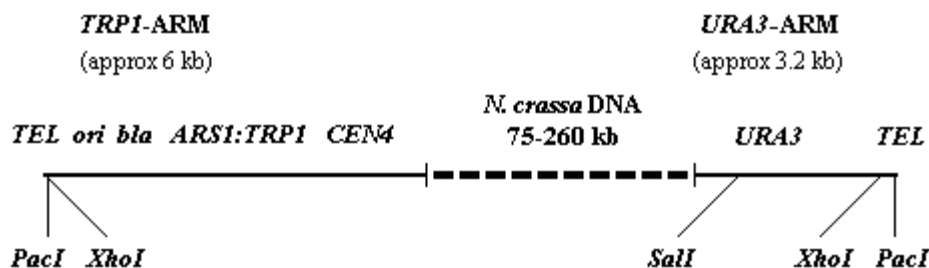


Figure 1. A schematic of a typical clone from the *N. crassa* YAC library is shown. The markers shown are described in the text.

The following wells in the library contain cultures that will not grow or grow slowly upon inoculation: 1:6A 1:11H 2:8F 2:11E 5:9G 6:5F 6:11E 10:1F 10:6A 10:9B 11:1H 11:2G 11:6A 12:7D 12:9C 14:3D 14:3E 14:9E 17:3G 17:6F 17:9H 17:10H 18:3E 19:7G 19:9F 20:2D 20:5E 20:7H 21:4F 22:3D 22:4A 22:9C 23:11D 24:4A 24:4B 24:4C. There are no clones in columns 5 to 12 on dish 24. Additionally, well 21:3E was contaminated with cells from well 21:3D. At least 5% of the clones besides those listed did not grow or grew poorly upon inoculation from a third generation copy of the YAC library distributed by the Fungal Genetics Stock Center (C. Yamashiro, unpublished observations).

Chromosome Walking

1. A protocol is outlined for hybridization screening of the YAC library. Positive YAC clones obtained from screening the library must be verified by Southern analysis. Mixed cultures within an individual well can occur; therefore, streak out cultures for single colonies on selective plates. Prepare genomic DNA from one or more individual colonies, and perform genomic Southern blot hybridizations with these DNAs as outlined below.

2. Several of the YAC clones obtained from a given library screening should be physically mapped with restriction endonucleases and the resultant maps compared. They should fit into a single set of contiguously overlapping sequences referred to as a "contig". Construction of a contig map assures the continuity of individual clones and avoids possible confusion resulting from multiple insert clones, whose restriction maps will not fit into the consensus maps generated from the bona fide majority of clones.

3. DNA fragments containing the terminal ends of the inserts of YAC clones which map to the periphery of the contig can be subcloned into *E. coli* or amplified directly (as described below). These fragments can then be used as hybridization probes for chromosome walking, or RFLP mapping.

Hybridization screening of the YAC library. (B. Brownstein et al. 1989. *Science* 244: 1348-1351, and B. Brownstein personal communication with modifications)

I. Preparation of a YAC library filter set.

1. Positively charged nylon membranes (MSI Inc, Westboro, MA) are layered onto the surface of large petri plates containing SD medium with casamino acids and adenine, "SD + CAA + A" (a supplemented minimal medium). Cultures from microtiter wells are replicated onto marked membranes using a multi-pin replicator. Each microtiter plate should be replicated onto a single membrane. Growth will take 3 to 4 days at 30°C. Fast growing colonies can be allowed to grow quite large (0.5 mm or larger) to allow the slow growing colonies to be well established (other protocols prefer smaller colonies so that cells are actively growing but lysis of old colonies works well here).

2. Prepare spheroplasting buffer CES/DTT/YLE (CDY). For 40 ml: Add 80 mg yeast lytic enzyme, 70,000 units/gram (ICN Biochemicals Irvine, CA) and 300 ml of 2 M dithiothreitol (DTT) to 40 ml CES solution (1 M sorbitol 100 mM Sodium citrate pH 7.0, 50 mM EDTA pH 8.0).

3. Place Whatman 3 MM or equivalent in lid of 150 mm Petri dish. Add approximately 7.0 ml CDY to Whatman and remove any bubbles (only a small amount of excess liquid should be present in lid). Carefully place a membrane containing YAC colonies onto the CDY saturated paper. Carefully remove bubbles from under the membrane.

4. Place the bottom of the Petri dish over the membrane and seal dish well with parafilm. Incubate 2 days 30-32°C.

5. Colony lysis and DNA denaturation is done by placing the membranes sequentially on Whatman 3MM paper saturated with the following solutions at room temperature for the times shown:

- a. 10% SDS/100 ug/ml proteinase K for 10 min.
- b. 0.5 N NaOH for 10 min.
- c. Three washes of 200 mM Tris pH 7.5/2X SSC for 5 min each.

6. Air dry the membranes on 3 MM paper for 2 h minimum (can go overnight). Bake for 2 h in vacuum oven at 80°C. Wrap membranes in plastic wrap (do not stack membranes, arrange them side to side). Store at -20°C.

II. Hybridization.

1. Prehybridize filters overnight at 65°C in a prehybridization (prehyb) solution of 1 M NaCl, 1% SDS and 10% dextran sulfate (Oncor Inc. Gaithersburg, MD). Four sets of six filters were each prehybridized with 50 ml of prehyb solution. Filters were incubated in plastic "Tupperware" containers with the lids on (hybridization bags could also be used). Discard the prehyb solution if cell debris is present (lots of cell debris will be present after the first prehybridization of the filter set).

2. Hybridization: Use 1 to 5 x 10⁵ cpm of 32P-labeled probe/ml of prehyb solution. Boil carrier salmon or herring sperm ssDNA (100 ug/ml in final hybridization) with the probe for 10 min to reduce background caused by non-specific hybridization. Add labeled probe with carrier DNA to hyb cocktail in fresh Tupperware containers then add filters from prehybridization. Four sets of six filters were each hybridized with 50 ml of cocktail. Hybridize at 65°C (use from 42°C to 65°C depending on stringency) overnight with shaking. To save money, we hybridized two sets of six filters one night, then reused the hybridization cocktail for the remaining filters the following day.

Note: YAC vector pYAC4 contains the entire pBR322 DNA sequence. Do not use probes with this sequence, otherwise your probe will hybridize to all the YAC clones. Alternatively, a 50-fold molar excess of cold competitor pBR322 DNA can be added to the hybridization reaction to reduce background. Similarly, if labeled cosmid DNA is used as a probe use a 50-fold molar excess of cold cosmid vector (i.e., lacking insert DNA) to decrease background (vector-vector) hybridization.

3. Remove membranes for washes on a gyratory shaker. Wash times should be adjusted to yield background signal acceptable to you.

- a. Wash in 2 changes of 2X SSC/ 1% SDS for 30 min each at 65°C.
- b. Wash in 2 changes of 0.5X SSC/ 1% SDS for 30 min at 65°C. For higher stringency use 0.1X SSC/ 1% SDS for 1 h at 75°C.

4. Wrap membranes in plastic wrap and expose to film.

5. Rehybridization: We did not strip our library filters. Instead, we prepared two filter sets and each set was allowed to decay approximately 3 months prior to rehybridization. If filters are to be stripped do not allow them to dry out at any time (store wet at -20°C).

Selective medium (SD + CAA +A) 0.67% Bacto-yeast nitrogen base without amino acids, 2% dextrose, 0.5% casamino acids, 15 mg adenine/liter. For solid medium add 2% agar.

YAC restriction mapping and contig building.

Restriction maps of YAC clones can be generated quickly by partial digestion and indirect end labeling (M. Centola and J. Carbon. 1994. *Mol. Cell. Biol.* 14:1510-1519) as follows:

Chromosomal DNAs of the YAC clones imbedded in agarose "DNA plugs" are prepared. DNA plugs are partially digested with a restriction endonuclease which recognizes an 8 bp DNA sequence. Partial digestion reactions are fractionated on CHEF gels and blotted to membranes. Membranes are then hybridized with probes specific for one of the YAC vector arms. A series of positive hybridization bands of increasing size will be imaged upon autoradiography of the membrane in each lane containing a partial digestion reaction. The size of each band is calculated relative to size standards run on the fractionation gel. The difference in size between sequential bands is equal to the distance between restriction sites in a given clone. Resultant restriction maps can be verified by rehybridizing stripped blots with a probe specific for the opposite arm of the YAC clone. Restriction maps of individual clones can then be overlapped and contigs can be generated.

I. Preparation of chromosomal DNA plugs of YAC clones (Kuspa et al., 1989. *Proc. Natl. Acad. Sci. USA* 86:8917-8921).

A. Small scale preparation:

1. Grow 3 ml culture of yeast strain in SD+CAA+A medium at 30°C to stationary phase (2-3 days).
2. Pellet 1.5 ml of cell culture in a microfuge tube and wash pellet with 1 ml 0.5 M EDTA pH 8.0.
3. Remove EDTA solution, respin cells and remove ALL remaining EDTA solution with a pipetman.
4. Gently resuspend cells in 10.8 µl spheroplasting solution.
5. Incubate at room temperature (RT) for 5 min.
6. Add 25.2 µl of molten 1.2% LMP agarose in SCE cooled to 45°C (final concentration of agarose is 0.6%).
7. Pour into plug molds. Solidify at RT for approximately 20 min.
8. Move plugs into microfuge tube. Add overlay solution to fill tube. Incubate overnight at 37°C to produce spheroplasts.

9. Remove overlay solution carefully.
10. Add 1 ml of NDS; incubate for 2 days at 50°C.

B. Large Scale Preparation

1. Grow 50 ml culture of yeast strain in a 125 ml flask, with shaking, to stationary phase (2-3 days) at 30°C in SD+CAA+A medium; cultures will turn pink.
2. Pellet cells 5 min at approximately 1000 x g.
3. Wash with 100 ml 50 mM EDTA pH 8.0 (EDTA solution).
4. Remove EDTA solution, respin cells and with a pipetman remove ALL remaining EDTA solution.
5. Gently resuspend pellet in 400 µl spheroplast solution.
6. Incubate at RT 5 min: briefly warm cells to 45°C (10-20 sec).
7. Add 800 µl 1.2% LMP agarose in SCE. Cool molten agarose to 45-48°C prior to addition. Mix gently and completely. Quickly pour into plug molds. Solidify at RT for approximately 20 min.
8. Move plugs into 15 ml Corex tube add 12 ml overlay solution. Incubate overnight at 37°C to produce spheroplasts.
9. Remove overlay solution carefully.
10. Add 10 ml of NDS and incubate overnight (for 1-2 days) at 50°C. Plugs can be stored in NDS at 4°C.

Solutions with volumes for large scale preps are given below. Final volumes can be scaled up or down.

Spheroplast solution

1 ml SCE
25 µl 2-mercaptoethanol
1 mg Zymolyase (100T) (ICN Biochemical, CA)

1.2% LMP agarose
0.12 g LMP agarose
10 ml SCE (melt agarose)

Overlay solution

46 ml 0.5 M EDTA

0.25 ml 2M Tris pH 8.0
3.75 ml 2-mercaptoethanol

NDS
1 mg/ml Proteinase K
1% N-lauroyl sarcosine in 0.5 M EDTA pH 9.0

SCE
1 M sorbitol
0.1 M sodium citrate, pH 5.8
10 mM EDTA (filter sterilize)

Plug molds: Use 1 cc syringe barrels with the ends sealed with parafilm or slide holders that come with common microscope slides. Tape up sides of slide holder and pour in molten agarose to form plugs that fit nicely into gel wells.

II. Partial restriction enzyme digestion of YAC DNA plugs.

Plugs must be extensively washed prior to digestion. Wash a given volume of plugs 5 times with 25 volumes of 10 mM Tris pH 8.0, 50 mM EDTA at RT for 2 hr/wash. Plugs can be washed and then stored in wash solution at 4°C for months prior to digestion.

1. Add the equivalent of approximately 20 μ l volume of solid chromosomal DNA plug/tube to two microfuge tubes (approximately 0.2-1 μ g of DNA). Two digestion reactions will be done on each clone.
2. To each tube add 100 μ l 0.5 mM phenylmethylsulfonyl fluoride (PMSF; USB, Cleveland, OH) in TE buffer. Incubate reactions at 37°C for 30 min. PMSF is highly unstable in aqueous solution (half life = 30 min); prepare a fresh 100 mM stock in 100% ethanol and dilute to 0.5 mM immediately before each use.
3. Remove PMSF wash solution and add 1 ml of TE buffer to plug. Incubate at 37°C for 30 min and repeat wash once. Remove final wash solution completely.
4. Add 1/10 volume of restriction enzyme buffer (10 x concentrate) to plug and melt mixture at 68°C for 5 min. Plugs should be completely molten.
5. Cool plugs to 37°C and add 0.2 and 1.0 units of pre-warmed (37°C) restriction enzyme to each of the two reactions, respectively. The plug should still be molten. Mix gently with a pipet tip. Incubate at 37°C for 30 min.
6. Add 1 μ l of 0.5 M EDTA (pH 8.0) and remelt plugs at 68°C for 5 min and the molten reactions loaded onto an agarose gel (see below).

Note: Complete restriction enzyme digests of the DNA within agarose plugs can be done as described above if a large amount of enzyme is used (we typically use 12-30 units of enzyme/20 μ l of plug) and incubate reactions at 37°C for 3 hr or more.

III. Using CHEF gel analysis to resolve YAC clones.

Partial digests are fractionated on a 1% agarose gel in 0.5 x TBE (1 x TBE = 90 mM Tris-borate, 2 mM EDTA). Gels (10 x 10 cm) are subjected to electrophoresis at 6.0 V/cm for 15 h in a CHEF gel apparatus (BioRad, Richmond CA), using a 0.2 to 13 sec pulse ramp, with a 120° pulse angle, at 14°C. Lambda concatemers and a 5 kb ladder (pBR328 concatemer, BioRad) size standards are also loaded on the gel. Stain gel with ethidium bromide, and photograph with a ruler atop the gel so migration distances of standards can be determined. These conditions will resolve DNA fragments ranging from 5-250 kb and are thus suitable for both partial digest mapping and for analysis of complete digestion products by Southern hybridization.

Some YAC clones are >250 kb. Resolution of uncut YAC clones can be done using the following conditions: 1% agarose gels prepared and subjected to electrophoresis in 0.5 x TBE buffer at 6.0 V/cm, using a 6-12 sec ramped switch time, at 12°C for 16-20 hr. These conditions will resolve clones up to 400 kb in size. Plugs can be placed into the wells of a CHEF gel by balancing the plug on a glass coverslip and pushing the plug into the well with a pipet tip. Alternatively, as suggested by Virginia Pollard, plugs can be placed onto the teeth of the well comb and the agarose gel cast with the comb in place. One needs to check for any dislodged plugs which can be repositioned using a Pasteur pipet.

IV. Southern Hybridization.

Southern blotting of high molecular weight DNA requires extensive nicking of the DNA prior to transfer, and a large volume of transfer buffer to be used during capillary transfer.

1. Nick high molecular weight DNA within the ethidium bromide stained gel using a commercial UV cross linker. Use settings recommended by the vendor; BioRad crosslinkers have a program for nicking pulse field gels that works nicely. Alternatively, place the gel into 0.2 N HCl for 30 min at RT with gentle agitation.
2. Denature gel in 0.5 M NaOH, 1.5 M NaCl for 30 min.
3. Neutralize in 1 M Tris pH 7.5, 1.5 M NaCl for 30 min.
4. Transfer to a positively charged nylon membrane using standard capillary blotting procedures with the following modifications:
 - a. Use a very large wick of 3MM paper, approximately 35 cm long x 30 cm wide, under the gel.
 - b. Use 2 liters of 10x SSC for each transfer. Use a large volume of absorbent material above the gel. Use two glass plates (approximately 30 cm x 20 cm) as weights above the absorbent, excessive weight will compress the gel and prematurely stop the transfer.
5. Standard Southern hybridization conditions can be used.
6. Hybridization probes specific for the left and right arms of the pYAC4 vector can be obtained from the following sources:

- a. Gel-purified 1.7 kb BamHI-PvuII fragment from pBR322 hybridizes specifically to the URA3-encoding YAC vector arm (URA3-arm).
- b. Intact pBluescript (Stratagene, La Jolla, CA) hybridizes specifically to the TRP1- encoding YAC vector arm (TRP1-arm).

Isolation of terminal restriction fragments from cloned DNA inserts in YAC clones.

Terminal restriction fragments of the insert DNA ("terminal fragments") from the TRP1-arm can be directly subcloned into *E. coli* by plasmid rescue and terminal fragments from the URA3-arm can be subcloned into *E. coli*. Alternatively, terminal restriction fragments from either end of YAC insert DNA can be cloned by ligation-mediated PCR. Although large regions of DNA (up to 30 kb) can be subcloned by plasmid rescue, the procedures are involved and at times frustrating. We therefore recommend the use of ligation-mediated PCR for obtaining small probes for chromosomal walking.

To isolate terminal end fragments of YAC inserts by plasmid rescue: Digest total genomic DNA from a given YAC clone with XhoI, or PacI (a double digest with XhoI and SalI should also work). These enzymes cut the vector arm once on the telomere proximal side of the *E. coli* bla and ori. These enzymes should also cut somewhere in the insert DNA. The digestion reaction is then diluted and ligated to circularize the DNA fragment containing bla, ori and the terminal restriction fragment of the insert DNA. A portion of the ligation reaction is electroporated into *E. coli* and AmpR colonies selected and subsequently screened for the rescued plasmid. Very large (>25 kb) plasmids do not rescue well into *E. coli*; therefore, determine the expected size of the rescue plasmid for a given enzyme before doing a rescue by preparing a Southern blot of CHEF gel fractionated genomic DNA cut with PacI, XhoI or XhoI/SalI from a given YAC clone and probing with a TRP1-arm specific probe. This will yield the approximate size of the terminal restriction fragment to be rescued. This information should be used to ensure that the rescued plasmid is of the correct size.

Similarly for subcloning terminal fragments from the URA3-arm, digest total genomic DNA from a given YAC clone with PacI, XhoI or XhoI/SalI and shotgun clone the entire genome into an *E. coli* vector (we used pBluescript). The yeast URA3 gene heterologously rescues the Ura phenotype of *E. coli* strain DB6656 lacZ 624(Am) trp-49(Am) pyrF79::Mu rpsL179 hsdR27 (Bach, et al. 1979 Proc. Natl. Acad. Sci. USA 76:386). The ligation reaction is introduced into this strain by electroporation and Ura⁺ transformants selected.

I. Preparation of miniprep DNA from liquid cultures of YAC clones (L. Clarke and M. Baum, personal communication).

1. Inoculate 10 ml of SD + CAA + A with yeast cells late in the day in a 125 ml or larger flask. Incubate overnight with shaking at 32°C until culture reaches an OD₆₆₀ of 0.6-0.7 (higher concentration cultures yield DNA resistant to restriction enzyme digestion). Cell cultures can be transferred to 50 ml disposable corning tubes upon reaching the proper OD, pelleted and stored at -20°C for several days before beginning DNA preps.

2. Wash with 2 ml 1 M sorbitol .

3. Loosen cell pellets before addition of SCE by brief vortexing, then resuspend cells gently in 2 ml SCE containing 0.2 mg Zymolyase 100T and 5 ul 2- mercaptoethanol. Incubate at 37oC for 1 h (never vortex the cells once Zymolyase is present).

4. Centrifuge cells at 800 x g for 4 min, and wash pellet with 2 ml SCE.

5. Remove wash, and resuspend cells in 1 ml NE solution (0.15 M NaCl, 0.1 M EDTA pH 8.0) containing 1% SDS. Add 5 ul 10 mg/ml DNase free RNase, incubate 15 min at 37oC.

6. Add 5 ul proteinase K (20 mg/ml) and incubate another 15 min at 37o

C. 7. Extract with 1.5 ml phenol/CHCl₃/isoamyl alcohol (25:24:1). Invert gently 10 times. Lay tubes on their sides for 5 min. Spin 10 min at 1600 g in clinical centrifuge. Remove aqueous layer with P1000 pipetman using pipet tips with cut off ends (cut off ends so orifice is approximately 0.3 cm in width). Interface is prone to mixing with aqueous phase so use caution. Repeat extraction once.

8. To aqueous layer add 0.12 ml of 7.5 M NH₄OAc and 1.2 ml isopropanol, invert to mix. Spin at 8500 rpm for 10 min.

9. Wash pellet with 5 ml of 70% ethanol. Air dry pellet, and resuspend in 50 ul TE buffer.

20 ml SCE: 10 ml 2 M sorbitol , 2 ml 1 M Sodium citrate pH 5.8, 2.4 ml 0.5 M EDTA pH 8.0, 5.6 ml H₂O.

Note: Be aware of the shear forces generated during the resuspension steps. Be gentle but thorough. Use cut off pipet tips when aliquoting DNA.

II. Plasmid rescue of the terminal restriction fragment from the TRP1-arm (L. Clarke and M. Baum, personal communication with modifications).

1. Digest 5 ul of miniprep DNA in a total reaction volume of 50 ul with 20-30 units of restriction enzyme for 3 h at 37oC. Miniprep yields will vary - 5 ul will yield about 2-5 ug. (save an aliquot of the digest to run on a pulse-field gel to confirm complete digestion and expected size as described above).

2. Add 40 ul H₂O + 10 ul 3 M NaOAc + 40 ul phenol/CHCl₃/isoamyl alcohol (25:24:1), "PCI" and mix.

3. Microfuge for 5 min to separate phases and transfer aqueous phase to a new tube.

4. Precipitate with 2 volumes of ethanol at -20oC for 1.5 h.

5. Spin in microfuge at full speed for 20 min.

6. Wash pellet with 0.5 ml 70% ethanol.

7. Resuspend pellet in 650 ul TE buffer and reprecipitate with 100 ul 7.5 M NH₄OAc + one volume of isopropanol (second precipitation helps to further clean up the reaction).
8. Pellet and wash as above.
9. Resuspend pellet in 1 ml of 1 x T4 DNA ligase buffer and add 2 units of T4 DNA ligase, incubate at room temperature overnight.
10. Split the reaction into two 500 ul aliquots; to each add 50 ul 3M NaOAc + two volumes of ethanol.
11. Store one aliquot at -20°C (back-up reaction) and the other precipitate at -70°C for 20 min.
12. Pellet and wash as before.
13. Resuspend pellet in 25 ul TE buffer.
14. Introduce a 1 ul aliquot of the reaction mix into *E. coli* by electroporation and select for Amp^R colonies. Note: Plasmid rescues containing repeats will tend to delete if ampicillin concentration in the selective media is >35 mg/l. Standard transformations can be used; however, it is more difficult to get rescues with lower transformation efficiencies.

III. Subcloning of terminal restriction fragment from the URA3 arm.

1. Digest 1 ug of miniprep yeast DNA with 20-30 units of restriction enzyme in a 50 ul final volume for 3 hr at 37°C as described above.
2. Digest 1 ug of pBluescript (or appropriate vector) with 20 units of enzyme as in step 1.
3. Add 0.5 unit of calf intestinal alkaline phosphatase (Boehringer Mannheim, Indianapolis, IN) to the pBluescript reaction, incubate an additional 20 min at 37°C.
4. Extract each reaction with an equal volume of PCI.
5. Precipitate DNAs with 1/10 volume of 7.5 M NH₄OAc + 1 volume of isopropanol.
6. Resuspend each reaction in 10 ul TE buffer, and check the DNA concentration by gel electrophoresis. (Alternatively, if you feel lucky assume your yield and proceed).
7. Ligate 100 ng of digested genomic DNA and 100 ng of prepared vector in a 20 ul reaction volume with 1-2 units of T4 DNA ligase overnight at 16°C.
8. Introduce 1 ul of the ligation reaction into *E. coli* strain DB6656 by electroporation and select for Ura⁺ transformants on *E. coli* minimal medium supplemented with 50 ug/ml tryptophan, and 50 ug/ml ampicillin (T. Maniatis et al. 1982. *Molecular cloning*, N.Y.). Colonies will appear after 2-3 days.

IV. Ligation-mediated PCR.

Ligation-mediated PCR (LM-PCR; Mueller and Wold 1989. Science 246: 780-786) is a procedure to isolate DNA fragments for which the sequence of the DNA is unknown. This technique has been adapted for generating DNA fragments corresponding to the end of inserts within YAC clones (Kere et al. 1992 Genomics 14: 241-248). Below we describe the LM-PCR technique used for generating end fragments of inserts within YAC clones with several minor modifications.

1. Six oligonucleotides are needed (5' to 3'):

L primer further from the cloning site (oligo L1)
CACCCGTTCTCGGAGCACTGTCCGACCGC

R Primer further from the cloning site (oligo R1)
ATATAGGCGCCAGCAACCGCACCTGTGGCG

L primer nearer the cloning site (oligo L2)
TCTCGGTAGCCAAGTTGGTTTAAGG

R primer nearer the cloning site (oligo R2)
GTCGAACGCCCGATCTCAAGATTAC

Linker long strand and linker primer
GCGGTGACCCGGGAGATCTGAATTC

Linker short strand
GAATTCAGATC

2. Digest 1 ug of yeast DNA embedded in agarose plugs (approximately 20 ul of plug) to completion as described above using restriction enzymes leaving blunt ends, such as *RsaI*, *AluI*, *PvuII*, *EcoRV*, or *ScaI* (usually set up 2-4 different digests for each clone). One can also digest 1 ug of total genomic DNA obtained from the miniprep procedure (above) in a 20 ul reaction volume.
3. Use 5 ul of the restriction digest and ligate the linker by adding ligation buffer, 25 pmol linker, and 1-2 units of T4 DNA ligase in a total reaction volume of 20 ul. Leave at room temperature for >1 h.
4. Set up the PCR using 2-5 ul ligation reaction, 20 pmol long vector-specific primer (L1 or R1) and 10 pmol of linker primer in a total volume of 50 ul. Cycling conditions are 94oC for 1 min, 65oC for 2 min, and 72oC for 2 min for 30-35 cycles.
5. Visualize products by subjecting 5-10 ul to electrophoresis on a 1.5% agarose gel.

6. To make probes from these fragments (they still have significant amount of vector sequence at this time), use PCR amplification of the purified first PCR product with 20 pmol of short vector-specific primer (L2 or R2) and 10 pmol of linker primer in a total volume of 50 μ l with the same cycle conditions used in step 4. Purifying the first PCR product is simple by first coring out a small amount of agarose containing the band of interest with a pasteur pipet and transferring to a microfuge tube. Add 100 μ l of TE buffer and boil for 5 min. Use 1 μ l for the second LM-PCR reaction. The PCR reaction can be done in the presence of 32 P-nucleotide for a radioactively labeled probe or in the presence of digoxigenin labeled dUTP (Boehringer Mannheim) for a non-radioactive probe (the latter probe can be used according to protocols provided by Boehringer Mannheim with the Genius kits and systems for non-radioactive detection of nucleic acids).

Transformation of *N. crassa* with a YAC clone.

I. Introduction of a *N. crassa* selectable marker into a YAC clone.

The pYAC4 vector does not encode a *N. crassa* selectable marker. Therefore, transformation of a given YAC clone is limited to clones that encompass selectable markers and must be carried out in appropriate mutant strains. Introduction of a clone from the library into *N. crassa* has been completed by transforming a qa-2, arom-9, inl strain of *N. crassa* (FGSC #3952) with total genomic DNA from YAC clone AB1380/YAC 12-10-H (M. Centola unpublished observations). This clone encompasses the qa-2+ gene. Standard spheroplast transformation procedures were used (Vollmer and Yanofsky 1986. Proc. Natl. Acad. Sci. USA 83:4869-4873) and 10 transformations using 1 μ g of DNA/150 μ l of competent cells yielded a single transformant.

To facilitate introduction of any YAC clone in the library into *N. crassa* a yeast integration vector carrying a dominant *N. crassa* selectable marker was constructed. The plasmid pLUSph (diagrammed in figure 2A) encodes the hygromycin phosphotransferase gene (hph) under the control of a modified cpc-1 regulatory region (Royer and Yamashiro 1992. Fungal Genetics Newsl. 39:76-79). This plasmid (pLUSph) also contains the yeast LYS2 gene, and an *E. coli* KanR gene and ori. This plasmid was constructed by inserting the 3.0 kb HindIII fragment from pMP6 into the unique HindIII site on pLUS (Hermanson et al 1991. Nucl. Acids Res. 18:4943-4948).

Integration of the pLUSph vector into the YAC clones

This plasmid can be site-direct integrated into the URA3-arm of the YAC clones by digesting pLUSph at the unique SalI site and introducing the linearized plasmid into the yeast host strain (diagrammed in figure 2B).

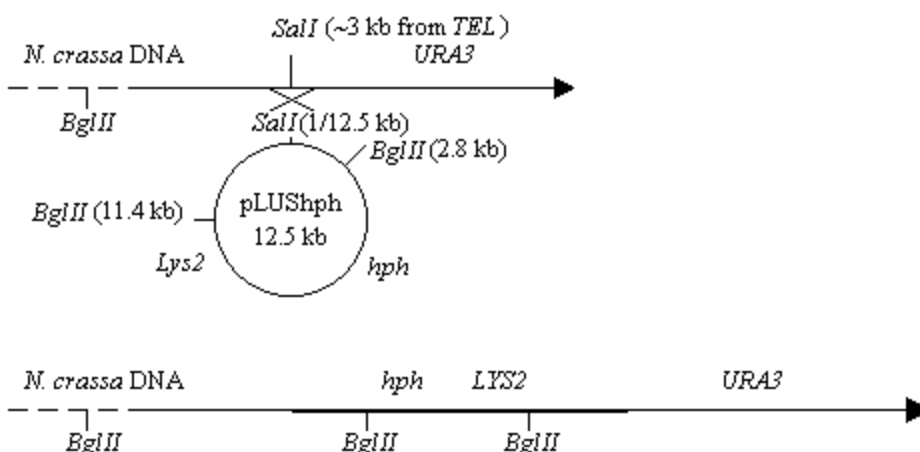


Figure 2. YAC integration vector pLUSph. A) A schematic of the YAC integration vector is shown. The region designated sup4-*URA3* is homologous to a portion of the *URA3*-arm of the YAC clones. As described in the text, this region can be utilized to mediate site-specific integration of this vector into the *URA3*-arm of any YAC clone. B) A schematic of the integration event is shown. The product of integration contains the yeast *LYS2* gene for phenotypic selection of integrated clones. A typical YAC integrant, shown at the bottom of the figure, also carries the *N. crassa* dominant selectable marker *hph* that can be used for selection of the YAC integrants upon introduction into *N. crassa*.

We have integrated this plasmid into several YAC clones by transforming competent yeast spheroplasts. However, lithium acetate transformations (Elble 1992. *BioTechniques* 13:18-20) are much easier to perform and transformation of LiOAc competent yeast cells with 1 ug of pLUSph linearized with *Sal*I should also be successful. The integrants are selected on dropout media lacking lysine and tryptophan (see below). The DNA from the colonies displaying the Lys⁺ phenotype was prepared in agarose plugs.

Dropout medium - SD + Adenine + essential amino acids (no lysine)

0.67% Bacto-yeast nitrogen base without amino acids

2% Dextrose

15 mg adenine/liter

2% agar.

Supplemented with the following amino acids. Final concentrations are shown in ug/ml.

aa (final conc.)

Arg (20)

Asp (100)

Met (20)

Glu (100)

Phe (50)

Val (150)
His (20)
Ser (375)
Tyr (30)
Leu (60)
Thr (200)

II. Verification of Integration by Southern hybridization and Restriction Enzyme Analysis.

A. Detection of an integrated copy of the hph gene within YAC clones by Southern analysis

1. Prepare DNA from the Lys⁺ yeast colonies in agarose plugs (see above for plug preparation protocol).
2. Fractionate the undigested DNA from the Lys⁺ transformants and parental host strains on a CHEF gel (use conditions described above for 5-250 kb resolution).
3. Transfer the DNA to a nitrocellulose membrane.
4. Hybridize the membrane with a radioactively-labeled DNA fragment containing the hph gene. We used the 1.2 kb BamHI/ClaI DNA restriction fragment from pCSN44 that encompasses the hph coding region (available from the FGSC). YAC clones in which pLUS_{hph} has integrated will hybridize to the probe.

B. Restriction Mapping of hph containing YAC clones.

1. Digest the DNAs from Lys⁺ YAC clones that hybridize to the DNA fragment containing the hph gene and the parental host strains with a restriction endonuclease which cuts inside pLUS_{hph}, but not inside pYAC4 (see "Restriction enzyme digestion of YAC DNA plugs" protocol shown above). We used the restriction endonuclease BglII.
2. Fractionate the digested DNAs on a CHEF gel (see "Using CHEF gel analysis to resolve YAC clones" protocol).
3. Transfer the DNA to a nitrocellulose membrane.
4. Probe the Southern blot with a radioactively-labeled DNA fragment which has homology with the YAC arm and pLUS_{hph}. We probed the Southern blot with the 1.7 kb BamHI/PvuII fragment from pBR322.
5. Expose the blot to film.
6. From the autoradiograph, determine the distance from the end of the telomere to the first restriction site within the insert of the parental YAC clone.
7. From the known restriction map of pLUS_{hph}, determine the size of the restriction fragments one would get if the plasmid were integrated at the unique SalI site. This is done by adding the distance from the end of the telomere to the SalI site on the YAC arm into which the plasmid was integrated to the distance between the SalI site on pLUS_{hph} and the first restriction site on

the SUP4 side of the plasmid. Conversely, one will see a band with a size corresponding to the distance from the restriction site within the insert to the SalI site on pLUSph to the next restriction site on the URA3 side of pLUSph.

Acknowledgements. We are grateful for the support and encouragement of John Carbon, Charles Yanofsky and David Perkins. We thank David Catcheside, Jack Kinsey, and Ed Cambareri for critically reading the manuscript.